

CD140a (PDGFRA) Monoclonal Antibody (APA5), PerCP-eFluor 710, eBioscience™

Product Details	
Size	100 µg
Species Reactivity	Mouse
Host/Isotype	Rat / IgG2a, kappa
Recommended Isotype Control	Rat IgG2a kappa Isotype Control (eBR2a), PerCP-eFluor 710, eBioscience™
Class	Monoclonal
Type	Antibody
Clone	APA5
Conjugate	PerCP-eFluor™ 710
Form	Liquid
Concentration	0.2 mg/mL
Purification	Affinity chromatography
Storage buffer	PBS, pH 7.2
Contains	0.09% sodium azide
Storage conditions	4° C, store in dark, DO NOT FREEZE!
RRID	AB_2815148

Applications	Tested Dilution	Publications
Flow Cytometry (Flow)	0.125 µg/test	5 Publications

Product Specific Information

Description: The APA5 monoclonal antibody reacts with the mouse CD140a molecule, the alpha chain of the platelet-derived growth factor receptor (PDGF receptor). PDGFRA is a receptor tyrosine kinase that forms dimers on the surface upon ligand binding and phosphorylates substrates. Dimers of PDGFR consist of either homodimers of alpha/alpha, beta/beta, or heterodimers of alpha/beta and serve as a substrate for its kinase activity. CD140a is expressed by embryonic tissues and mesenchymal-derived cells of the adult mouse tissues.

Applications Reported: This APA5 antibody has been reported for use in flow cytometric analysis.

Applications Tested: This APA5 antibody has been tested by flow cytometric analysis of NIH/3T3 cells. This may be used at less than or equal to 0.125 µg per test. A test is defined as the amount (µg) of an antibody that will stain a cell sample in a final volume of 100 µL. Cell number should be determined empirically but can range from 10⁵ to 10⁸ cells/test. It is recommended that the antibody be carefully titrated for optimal performance in the assay of interest.


PerCP-eFluor 710 emits at 710 nm and is excited with the blue laser (488 nm); it can be used in place of PerCP-Cyanine5.5. We recommend using a 710/50 bandpass filter, however, the 695/40 bandpass filter is an acceptable alternative. Please make sure that your instrument is capable of detecting this fluorochrome.

Light sensitivity: This tandem dye is sensitive to photo-induced oxidation. Please protect this vial and stained samples from light.

Fixation: Samples can be stored in IC Fixation Buffer (Product # 00-8222-49) (100 μ L of cell sample + 100 μ L of IC Fixation Buffer) or 1-step Fix/Lyse Solution (Product # 00-5333-57) for up to 3 days in the dark at 4°C with minimal impact on brightness and FRET efficiency/compensation. Some generalizations regarding fluorophore performance after fixation can be made, but clone specific performance should be determined empirically.

Excitation: 488 nm; Emission: 710 nm; Laser: Blue Laser

Product Images For CD140a (PDGFRA) Monoclonal Antibody (APA5), PerCP-eFluor 710, eBioscience™

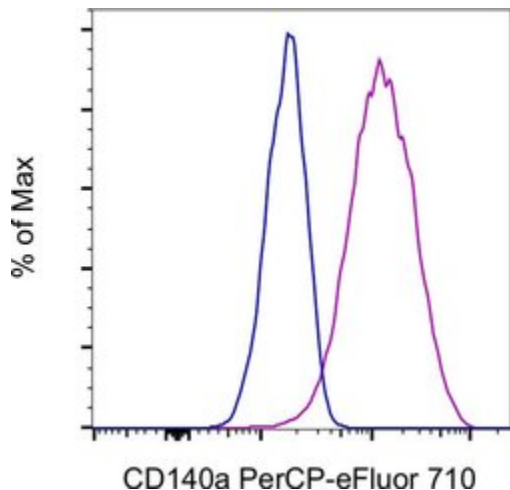

CD140a (PDGFRA) Antibody (46-1401-82) in Flow

NIH/3T3 cell line was stained with 0.06 μ g of Rat IgG2a kappa Isotype Control, PerCP-eFluor 710 (Product # 46-4321-82) (blue histogram) or 0.06 μ g of CD140a Monoclonal Antibody, PerCP-eFluor 710 (purple histogram). Total viable cells were used for analysis, as determined by Fixable Viability Dye eFluor 780 (Product # 65-0865-18).

[View more figures on thermofisher.com](https://www.thermofisher.com)

5 References

Flow Cytometry (5)

JCI insight

Cardiac fibroblast proliferation rates and collagen expression mature early and are unaltered with advancing age.

"Published figure using CD140a (PDGFRA) monoclonal antibody (Product # 46-1401-82) in Flow Cytometry"

Authors: Wu R, Ma F, Tosevska A, Farrell C, Pellegrini M, Deb A

Species

Not Applicable

Dilution

Not Cited

Year

2020

JCI insight

Loss of Fas signaling in fibroblasts impairs homeostatic fibrosis resolution and promotes persistent pulmonary fibrosis.

"Published figure using CD140a (PDGFRA) monoclonal antibody (Product # 46-1401-82) in Flow Cytometry"

Authors: Redente EF, Chakraborty S, Sajuthi S, Black BP, Edelman BL, Seibold MA, Riches DW

Species

Not Applicable

Dilution

Not Cited

Year

2020

[View more Flow references on thermofisher.com](#)

More applications with references on thermofisher.com

For Research Use Only. Not for use in diagnostic procedures. Not for resale without express authorization. Products are warranted to operate or perform substantially in conformance with published Product specifications in effect at the time of sale, as set forth in the Production documentation, specifications and/or accompanying package inserts ("Documentation"). No claim of suitability for use in applications regulated by FDA is made. The warranty provided herein is valid only when used by properly trained individuals. Unless otherwise stated in the Documentation, this warranty is limited to one year from date of shipment when the Product is subjected to normal, proper and intended usage. This warranty does not extend to anyone other than the Buyer. Any model or sample furnished to Buyer is merely illustrative of the general type and quality of goods and does not represent that any Product will conform to such model or sample. NO OTHER WARRANTIES, EXPRESS OR IMPLIED, ARE GRANTED INCLUDING WITHOUT LIMITATION, IMPLIED WARRANTIES OF MERCHANTABILITY, FITNESS FOR ANY PARTICULAR PURPOSE, OR NON INFRINGEMENT. BUYER'S EXCLUSIVE REMEDY FOR NON-CONFORMING PRODUCTS DURING THE WARRANTY PERIOD IS LIMITED TO REPAIR, REPLACEMENT OF OR REFUND FOR THE NON-CONFORMING PRODUCT(S) AT SELLER'S SOLE OPTION. THERE IS NO OBLIGATION TO REPAIR, REPLACE OR REFUND FOR PRODUCTS AS THE RESULT OF (I) ACCIDENT, DISASTER OR EVENT OF FORCE MAJEURE, (II) MISUSE, FAULT OR NEGLIGENCE OF OR BY BUYER, (III) USE OF THE PRODUCTS IN A MANNER FOR WHICH THEY WERE NOT DESIGNED, OR (IV) IMPROPER STORAGE AND HANDLING OF THE PRODUCTS. Unless otherwise expressly stated on the Product or in the documentation accompanying the Product, the Product is intended for research only and is not to be used for any other purpose, including without limitation, unauthorized commercial uses, in vitro diagnostic uses, ex vivo or in vivo therapeutic uses, or any type of consumption by or application to human or animals.